


Surpahs® W1 Body Weight Scale User Manual (V1.0)

The user manual can be downloaded through Surpahs® website at <http://www.surpahs.com>




Scan for shopping and support:

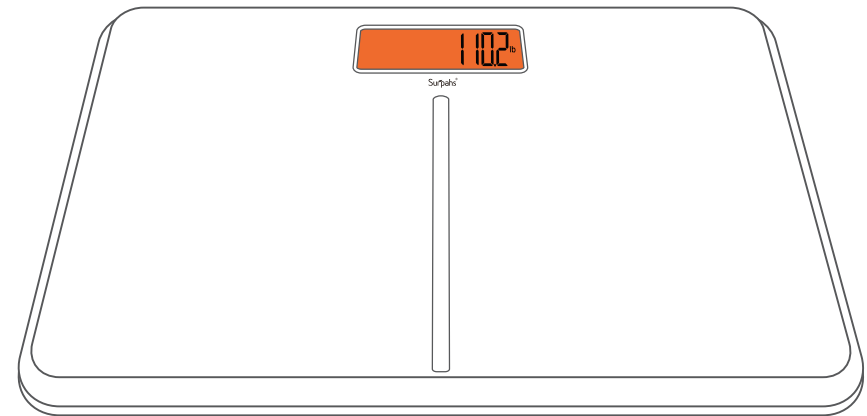

Surpahs®
www.surpahs.com

FEATURES & SPECIFICATIONS


Dimension	310 x 300 x 23.5mm (12.2" x 11.8" x 0.93")
Turn On Method	Sense-ON (Stand on the platform directly)
Sensors	4 high precision G-Sensors
Display	3.5" Large LCD with Orange Backlight V.A.=76mm (L) x 40mm (W)
Capacity	11 lb/5 kg - Max. 397 lb/180 kg
Division	0.1kg ; 0.2 lb ; 0.2 lb/1 lb (≥ 20 st)
Units	Kilogram (kg) / Stone (st:lb) / Pound (lb)
Power	6V (4xAAA Batteries)
Auto-Off Time	About 10 seconds of Non-Operation
Accuracy	0-50kg (0-110lb): ± 0.3 kg (0.66 lb) 50-100kg (110-220lb): ± 0.4 kg (0.88 lb) 100-150kg (220-330lb): ± 0.5 kg (1.10 lb) 150-180kg (330-396lb): ± 0.7 kg (1.54 lb)
Working Environment	Temperature: 0°C to 40°C / Humidity: $\leq 90\%$ RH
Storage Environment	Temperature: -20°C to 60°C / Humidity: 10%RH to 93% RH
Low battery indication	"Lo"
Overload Indication	"----"

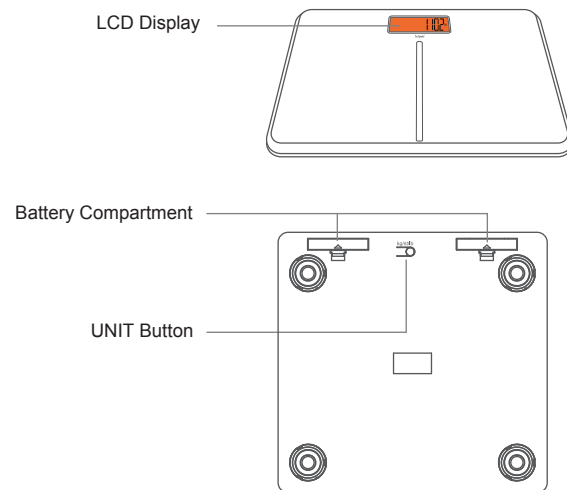

BATTERIES

- Uses 4 x AAA batteries
- Make sure the + and – sides of the batteries are aligned correctly.
- Replace batteries when the low battery indicator is shown.

When batteries are inserted, the screen will display "8888".

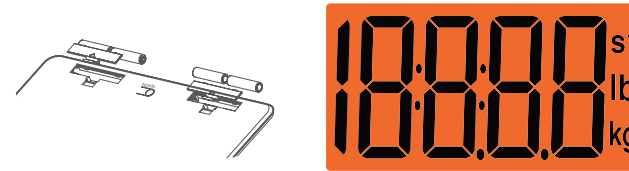

Wait until the screen display "0.0", and it is ready to use.



OPERATION INSTRUCTION

I. Initialization

NOTE: To ensure measuring accuracy, this initialization process must be done at anytime the scale is moved to another place, or batteries are changed.


- Press the platform center and remove your foot.
- "0.0" is displayed.



- The scale will switch off and now is ready to use.

II. Weighing

- Place the scale on hard flat even surface for accurate readings.
- Stand on the scale platform directly, it will turn on and measure your weight automatically

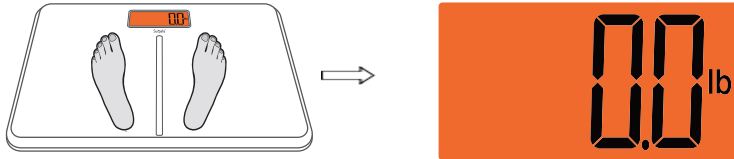

- Stand still and the weight data will be locked with twice blinking when the weight is stable. Read the result. Then the scale will turn off in about 10 seconds.



III. UNIT Selection

- To change measuring unit, turn the scale upside down and press the UNIT button until the desired unit is selected.
- You can choose from Kilogram (kg), Pounds (lbs) or Stones (st/lb)

IV. Turning Off

The screen will turn off in about 10 seconds.

TROUBLESHOOTING

Issues	Cause	Solution
	Overload. The scale will turn off in 4 seconds.	Stop using the scale for this measurement, or it will damage the G-Sensors.
	Low battery. The scale will turn off in 4 seconds.	Replace with new batteries

MAINTENANCE AND CARE

- Please use the scale on a flat, level and hard surface free from vibration and drafts, use it at room temperature.
- If the scale is not to be used for a long period, please remove the batteries,
- Clean the scale with a slightly damp cloth, do NOT immerse the scale body in water or use chemical/abrasive cleaning agents.
- Check the batteries first if you are having any trouble with your scale. This simple step can remedy most scale issues.
- Do not jump onto the scale platform immediately.
- Avoid storage in following locations:
 - Where there is water.
 - Where the device may be exposed to extreme temperatures, humidity, moisture, direct sunlight, dust, or salt air.
 - Where there is risk of shock or drop.
 - Where you store chemicals or full of corrosive gases.
 - Where in reach of the infants or children

WARRANTY

This product is warranted against defects in materials and workmanship for one year from the date purchase, when used in accordance with the instructions provided. This warranty does not cover damages or wear resulting from accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair.

C8KE INC. shall not be liable for loss of use or any other incidental, consequential or indirect costs, expenses or damages. There are no express warranties except as listed above. This warranty gives you specific legal rights, and you may have other rights which vary from state to state.

Should this product require service (or replacement at our option) while under warranty, please visit website www.surpahs.com and submit a support ticket.

SUPPORT

If you have any questions about this product, please visit www.surpahs.com, or scan the QR code, find the support link.

Surpahs®, C8KE INC.
DE 19809, USA
www.surpahs.com